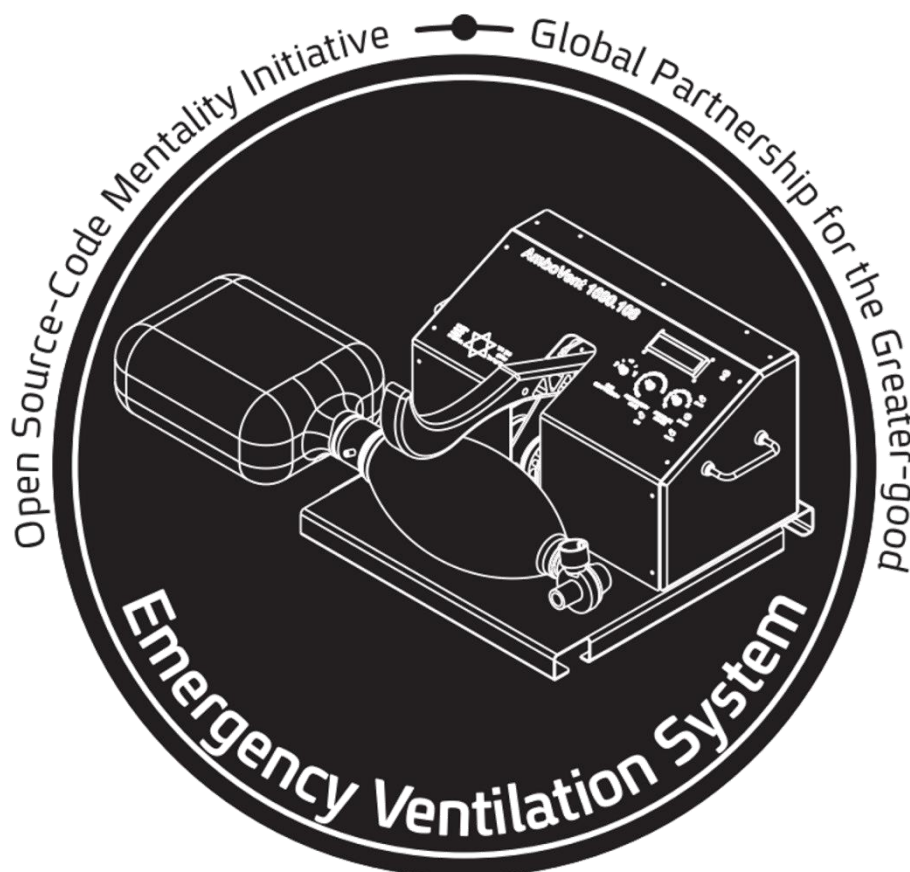


AmboVent
1690.108

Emergency ventilation alternative system

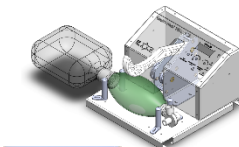
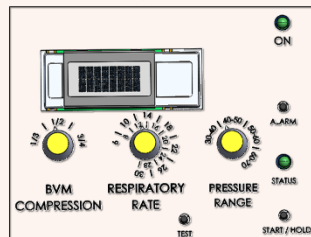
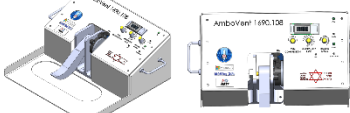
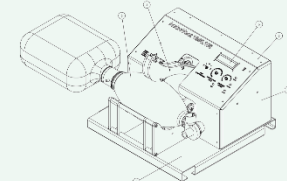
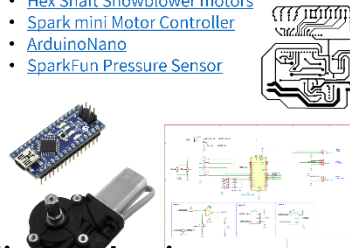


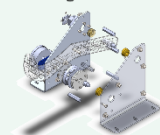



Open-Source Code Mentality Initiative
Global Partnership for the Greater-good

AmboVent 1690.108 for COVID-19 ventilation machines shortage challenge

1. Emergency Alternative Ventilation System: An open Source Initiative

Created by teams of Israeli volunteers: First Israel Robotics team, IAI Unit 108, Magen David Adom, Rafael, Microsoft R&D Israel, IDC Herzliya, and many mentors. Designed as an alternative automatic, controlled, ventilation system for adults, to be used only in emergency situations when no other approved ventilation machines are available.

<p>Specifications</p> <ul style="list-style-type: none"> • Volume control ventilation (VCV) machine • 3 tidal-volume respective to (33%, 50%, 75%) of the bag total volume • 11 ventilation rates ranging from 6 to 24 cycles per minute • 5 ventilation steps between 30-70 cmH2O, of not-to-exceed inspirium pressure • Electrical Source 110/220V • Batteries Backup (2 hours duration) • Visual and Vocal Alerts: Batt. On Low Batt Hose disc. Vent. rate fail. Extreme Pressure • Ventilation pressure continuous monitoring • IOT-based system, to monitor multiple patients 	<p>Emergency ventilation initiative coming out of Israel</p> <p>Who is behind it? Lead by the CTO & innovation leader of the Israeli Air-Force 108 Electronic Depot and backed by a large community of innovators behind him</p> <p>To include: 40 Professional Volunteers - Israel's national EMS, Physicians from leading Israeli hospitals such as Tel Aviv Sourasky and Hadassah JLM as well as other medical centers, Engineers, First Israel mentors and students, The Haifa Technological Center Rafael and Israel Aerospace Industries, IAF Unit 108, The garage program by Microsoft Israel and others..</p>	 <p>Analysis of Open Source COVID-19 Pandemic Ventilator Projects</p> <p>Rated by Robert Lee Reed as the leading solution in this category, among 40 other initiatives</p>
<p>UI-UX Specifications</p> 	<p>AmboVent 1690.108</p> <p>Automatic, Controlled Resuscitator Device</p>  <p>Emergency ventilation alternative system Global Partnership for the Greater good Leading Open Source Code Mentality Initiative</p>	<p>Documentation</p> <ul style="list-style-type: none"> • Files are available in an ANSI-metric standard • Mechanical model design source-code files are provided in a Parasolid format (.X_T) 
<p>Off the Shelf Components</p> <ul style="list-style-type: none"> • Hex Shaft Snowblower motors • Spark mini Motor Controller • ArduinoNano • SparkFun Pressure Sensor 	<p>Validation and Tests</p> <ul style="list-style-type: none"> • Supports Rapidly Manufactured Ventilator System specifications • Evaluated in Tel Aviv Sourasky hosp. • Tested in Flex labs • Calibrated in Rambam hospital  	<p>Getting Ready for Production</p> <ul style="list-style-type: none"> • Materials: Aluminum, Akulon-nylon6 • Metallic painting capacities • FDM based 3D printing machines • Factory assembling capabilities • CNC, Punch, Bend manufacturing Capabilities • Printing on Lexan  <p>Maj. Dr. David Alkaher Dalkaher@gmail.com</p> 

2. Medical evaluation

AmboVent 1690.108 is currently being tested by the Israeli Regulatory Authority. We work around the clock to complete both clinical performance and safety evaluation, including electromagnetic compatibility (EMC). So far, the testing shows the device can serve, in case there is no other regulatory-approved ventilator available, as an effective alternative ventilator for patients who require mechanical ventilation.

3. Technical requirements:

- Manual selection of tidal volume, by choosing % of squeeze from full (100%) bag squeeze. For example, choosing 60% means the device pushes out every cycle, 60% of its full (100%) capacity.



- b. Enable selection of 9 different respiratory rates, from 6-24 cycles per minute, increasing by steps of 2.
- c. Compatible with Ambu (and other manufacturers) commonly used self-inflating bags, ranging between 1100 to 1475 cc.
- d. Maximum inhalation pressure setting, ranging between 30 to 70 cmH₂O, increasing by steps of 10.
- e. Works with standard Positive End-Expiratory Pressure (PEEP) valves.
- f. Standard 110-220V Powerline feed + two-hours battery backup.
- g. In events of resistance during inhalation (abnormal rise in air pressure), the bag squeezing process stops and regains inhalation in the following cycle.
- h. Insulation of electrical and electronic parts to prevent the chance of sparking (in a pure oxygen-rich ventilation environment).

4. **Functionality:**

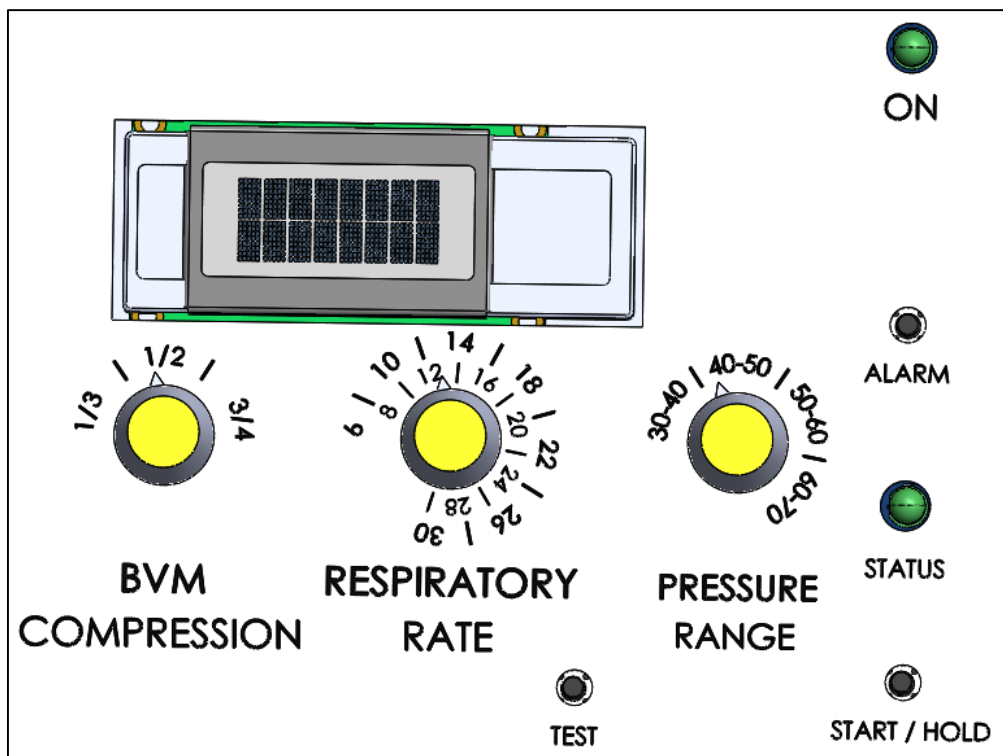
- a. Compliance with the use, connection, parts and integration of ventilation piping available in medical facilities.
- b. Compact and lightweight. Can be positioned with flexibility around the patient's bed, up to 1.5 meters away with no fear of increasing the dead space.
- c. Wide options of tidal volume selection, ranging from 30% to 100% of full squeezing capacity.
- d. Nine different respiratory rates, ranging between 6 to 24 cycles per minute.
- e. A predetermined I:E time ratio of 1:2. Can be changed in the program.
- f. Positive End-Expiratory Pressure (PEEP) control, using a standard PEEP valve.
- g. 110-220V power supply
- h. Two hours of continuous operation on backup batteries in case of external power supply failure.
- i. Simple, durable and intuitive structure and operation.
- j. Capable of choosing Peak Inspiratory Pressure (PIP) sensing threshold, ranging between 30-70 cmH₂O.
- k. Compatible with standard ventilation oxygen bags.
- l. Can be connected to hospital's clean air supply. Useful for ventilating fresh and clean air in crowded rooms.

5. **Automatic alerting upon malfunctioning**




- a. Electrical power supply failure (Audio alert (sound alarm) + visual indication by LED).
- b. Internal battery voltage drop (a continuous sound alarm that cannot be dismissed before reconnecting to an external power source + a visible led indication).
- c. Alerting in case someone tries to turn off the device during active ventilation (will cause a one-minute continuous sound alarm + a visible led indication).

- d. Pressure rise above the PIP threshold (one-time sound alarm. In case it continues (repeating events) – (continuous sound alarm with a prominent light indicator).
- e. Sudden, unexpected pressure drops (may indicate air tubes disconnect) (continuous sound alarm that stops if pressure is build back to normal within two ventilation cycles. Otherwise, the alarm continues until operator active intervention).
- f. Deviation from the user respiration rate setting (continuous sound alarm + a visible led indication).

User Interface and Operation



6. UI: Functionality

Display	Type	Remarks
BVM Compression	Potentiometer 	Pick values ranging from 30% - 100% of full bag squeeze
Respiratory Rate	Potentiometer	Pick values between 6 to 24 (increments of 2).
PIP Range	Potentiometer	Pick values between 30 to 70 cmH ₂ O (increments of 10).
Start / Hold	Button	Single short press turns the device on / 5 seconds continued press turns the device off.
Test	Button 	Pressing the button runs a one time cycle of 75% bag squeeze.
Alarm	Button	Turns off screen alerts and LED Status
Status	LED 	<p>LED Indications</p> <p>The device is on in Hold mode - Green</p> <p>Ventilation is on - Blinking green, synchronized with the respiration rate</p> <p>Malfunction - Blinking orange, synchronized with the respiration rate (error type will be listed on screen)</p> <p>Critical malfunction - Red</p> <p>Transitioning to battery mode (power disconnected) - blinking red in synchrony with the respiration rate</p>

On	LED	Indication of turning on the device - green in operating mode
-----------	-----	---

7. UI: Feedback for User

Display	Type	Display	Remarks
Low pressure	Number	Value of minimum pressure as measured by the system	Check the system
Ventilation pattern differs from operator settings	Number	Value of maximum pressure as measured by the system	Reenter value of preferred ventilation rate
Alert when working on battery power	Potentiometer	"Batt On" Indication after disconnecting power source	LED - blinking red
Alert when continually working on battery for 2 hours	Button	"Low Batt" Indication of reaching 2 hrs of continuous operation on batteries	LED - malfunction blinking orange
Alert when ventilation hose disconnects	Button	"Hose Disconnect" Indication pressure hose disconnected	LED - malfunction blinking orange

Emergency Ventilation Initiative Coming Out of Israel



Critical malfunction alert	Button	"Vent. Rate fail"	LED - malfunction red
Pressure reaching the PIP threshold	Number	Extreme Pressure	Alarm

Emergency Ventilation Initiative Coming Out of Israel



The following has volunteered & contributed to the success of this project

Email	Name	Last Name	Cell Phone Number	Job Title	Organization
Dalkaher@gmail.com	David	Alkaher	972542211442	CTO & Innovation Leader	-
shlomo@alkaher.com	Shlomo	Alkaher	972549990700	CEO	Alkaher Technologies
dreliram@gmail.com	Eitan	Eliram	972506563762	CEO	Interactive Innovation
kornblau@gmail.com	giora	Kornblau	972526606910	CEO	Terra ventures
zeev@shneorson.com	Zeev	Shneorson	972547232222	R&D manager	NovoCogito
yuvaleran@gmail.com	Yuval	Eran	972524764102	Paramedic, MDA project manager, team OnyxTronix executive mentor	MDA & FIRST
aryehkohn@gmail.com	Aryeh	Kohn	972545465905	CTO of medical technology and simulation center	Tel Aviv medical center- Ichilov
gmeiry@gmail.com	Dr. Gideon	Meiry	972545220015	Entrepreneur & Mentor in the Bio-Medical Field	CTO at Galil-Ofek Technological Incubator
yeru1970@gmail.com	Michael	Yerushalmi	972529452702	Mechanical engineer- Product Development mannager	
shtefandenis@gmail.com	Denis	Shtefan	972547459020	Mechanical Design Engineer	Israel Air Force - Base 108
nissim.sofer@gmail.com	Nissim	Soffer mor	972529469697	R&D technological consultant	108 IAF
ASHER.LEVIN@GMAIL.COM	Asher	Levin	972529434355	R&D team leader	IAF
palmedic@gmail.com	Guy	Shalev	972547775115	Manager of Microsoft innovation lab Israel	Microsoft
adargat1234@gmail.com	Adar	Gat	972528781584		
Lichi@lichi.co.il	Mikhal	Heffer	972523516070	Founder and CEO	Lichi Translations
kornblauy@gmail.com	Yair	Kornblau	972549414477	High School Student / Mechanical Design Lead	Orbit 1690
sagi8898@gmail.com	Sagi	Froylich	972504283905	Mechanical Engineer	OnyxTronix - Team 2231
noamalch@gmail.com	Noam	Alchanati	972544242174	First Robotics Competition Mentor	Steampunk 1577
GuyShneorson@gmail.com	Guy	Shneorson	972545676642	Mechanical&Material science Engineer	IDF R&D unit
arbib123@gmail.com	Jonathan	Arbib	972548090864	FIRST Alum, mentor in FRC team - Steampunk 1577,	FIRST - FRC Team - Steampunk 1577

Emergency Ventilation Initiative Coming Out of Israel



dortorg111@gmail.com	Dor	Turjeman	972526284757	Software team leader	IAF - 108
baruchba@gmail.com	Barouch	Bakalo	972545749259	Vice President Of Public Relations	Haifa Technological Center
orhirshfeld@gmail.com	Or	Hirshfeld	972545844555	Mentor	Haifa Tech Center
yossi.rabinovitz@gmail.com	Yossi	Rabinovitz	972546620317	CEO	Quay
kozirovsky@gmail.com	Sergey	Kozirovsky	972584240980	Mechanical designer	
alon@aran-rd.com	Alon	Nahshon	972546313706	VP R&D	ARAN Reaserch and Development
nimrodma99@gmail.com	Nimrod	Machlav	972509010197	Integration of software and electronics	FIRST FRC OnyxTronix 2231 at "Tichon Shoham" high school
osherbar123@gmail.com	Osher	Bar	972547605231		
tomkertesz10@gmail.com	Tom	Kertesz	972584231479	Mechanical Engineer	Orbit #1690
noa.singer.123@gmail.com	Noa	Singer	972587973145	Resource Coordinator	FIRST OnyxTronix
roee.baryadin@gmail.com	Roe	Bar-Yadin	972525413723	Mechanical designer	IDF R&D unit

The success of this project is also attributed to

[Israel's national EMS](#), [Tel Aviv Sourasky](#), [Hadassah JLM](#), [First Israel](#), [Haifa Technological Center](#), [Israel Aerospace Industries](#), [The garage program by Microsoft Israel](#), [IDC](#), [lichi](#), and others